

SYLFIRM X

Next Generation Dual Wave RF Micro Needling Platform

Unlike conventional RF microneedling systems, SYLFIRM X has both Continuous Wave and Pulsed Wave therefore capable of treating pigmented and vascular lesions such as melasma and diffuse redness, as well as skin rejuvenation and scarring. In addition, dual wave is subdivided into 8 different modes, which can be customised to individual patients, taking into account the indications and skin conditions to provide a delicate, precise and bespoke treatment.

INDICATIONS

FACE

- Pigmentation
- Melasma
- Post Inflammatory Hyperpigmentation (PIH)
- Vascular Lesions - Rosacea, Telangiectasia
- Skin Laxity
- Skin Rejuvenation
- Skin Resurfacing
- Facial Contouring
- Wrinkles & Fine Lines

BODY

- Stretch Marks
- Surgical Scars
- Atrophic Scars
- Sagging Skin



NOVUSMEDICAL

Clinically Proven Technology




NOVUSMEDICAL

Clinically Proven Technology


CONTINUOUS WAVE (CW) MODE

- Suitable for skin rejuvenation, tightening, lifting and scar treatment.
- CW mode creates wide, independent coagulation around each electrode.
- The coagulation area continues from deeper dermis to upper dermis with one penetration.
- Promising dermal remodeling and Skin Tightening.
- Clients can see the skin contraction immediately.

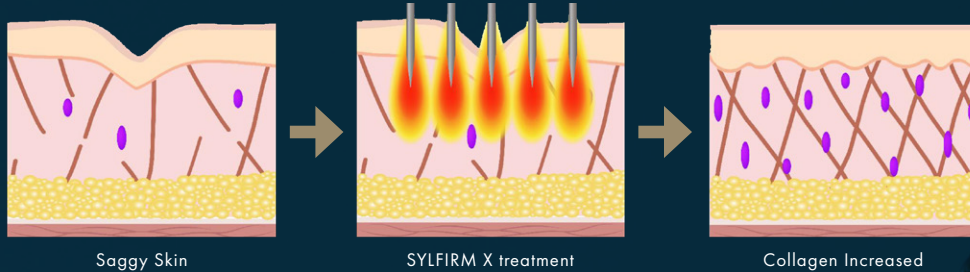

PULSED WAVE (PW) MODE

- Suitable for treatment of pigmented and vascular lesions such as melasma and diffuse redness.
- PW mode creates a soft coagulation selectively in basal layer & upper dermis.
- It coagulates B.M. and abnormal small blood vessels, while the surrounding tissue while the surrounding tissue is intact.
- Selective coagulation and dermal remodeling help to recover from hyper active melanogenesis as well as hyper angiogenesis.

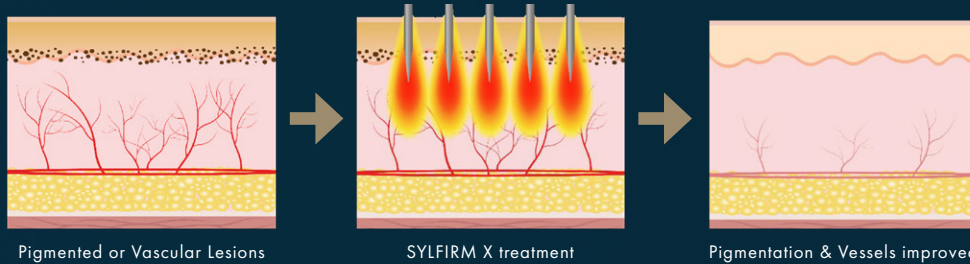

SYLFIRM X SUMMARY

- Continuous Wave and Pulsed Wave in one device
- 12 different treatment options
- Pigmentation, Melasma, Vascular Lesions, Rosacea, Telangiectasia (when used in Pulsed Wave mode)
- 25 gold plated non-insulated needles (thickness 300µm)

SPECIFICATIONS

SYSTEM TYPE	RADIO FREQUENCY MICRO-NEEDLING SYSTEM
ELECTRODE TYPE	BIPOLAR
FREQUENCY	2MHZ
MODES	CW1, 2, 3, 4 / PW1, 2, 3, 4
ENERGY INTENSITY	LEVEL 0-12
NEEDLE DEPTH	0.3-4.0 MM ADJUSTABLE *0.1MM STEP
NEEDLE SPEED	0.1-1.0 SEC ADJUSTABLE *0.1 SEC STEP
TIPTYPES	X TIP (5 X 5 PIN, 2MM SPACING) / XE TIP (5 X 5 PIN, 2MM SPACING) XA TIP (5 X 5 PIN, 1.5MM SPACING) / XW TIP (3 X 6 PIN, 1.5MM SPACING)
OPERATING METHOD	FOOT SWITCH
ELECTRICITY	AC 100 ~ 240V, 50/60HZ
DIMENSION / WEIGHT	440 W X 580 D X 1160 H CM / APPROX. 40KG
CERTIFICATES	FDA,CE, HEALTHCANADA, SFDA, KFDA

